

Seg-WAYTM ECuTR

Endoscopic Cubital Tunnel Release (ECuTR)

Seg-WAY is the first ECuTR system that allows the surgeon to safely view the entire ligament from above and below through controlled portals in the guide. It is also the only system in which the scope and the blade work independently of each other giving the surgeon a better field of view and a tactile feel. The system is designed so that the Ulnar nerve may be safely viewed and followed before the release is performed. The guide covers the blade on both sides of the ligament so as to protect any unwanted release.

Seg-WAYTM ECuTR

Features and Benefits

SAFETY FIRST

- Split Guide/Dual Scope Channels - 360° view allows the surgeon to view the ligament above and below to ensure there is no nerve
- Nerve Portal - allows the surgeon to rotate the scope to view and follow the nerve
- Covered Blade Channel - protects surrounding tissue and nerve while releasing the ligament

UNIPORTAL

- Smaller incision for faster recovery

VERSATILITY

- Rotating Guide - allows the surgeon to release with either hand or to move the release over

ASSURANCE

- Large field of view above and below to make sure no nerve is in the release

CONFIDENCE

- Longer guide allows the surgeon to release as much as necessary




Uniportal dual scope guide system



Individual long guide system



Independent antegrade knife operation with view port

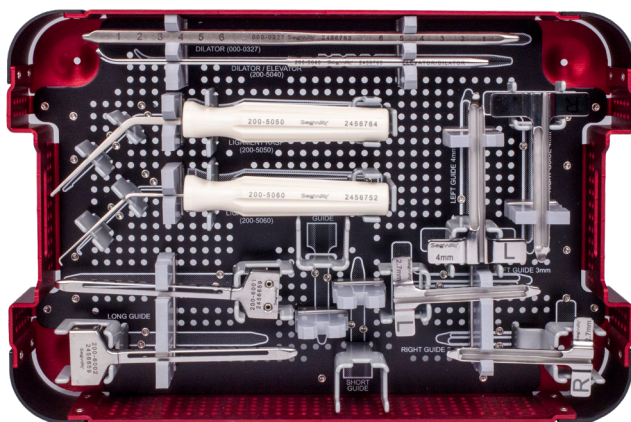

Description	Part #
Seg-WAY Instrument Set (Complete)	200-5560
Left Large SS Knife Guide	200-5012
Right Large SS Knife Guide	200-5022
Left Small SS Knife Guide	200-5014
Right Small SS Knife Guide	200-5024
Synovial Dilator / Elevator	200-5040
Rasp, Ligament	200-5050
Probe, Ligament	200-5060
Dilator	000-0327
Antegrade Ligament Knife (Single Use)	200-1006



40 General Warren Blvd, Suite 100 | Malvern, PA 19355 U.S.A.
610.989.8080 | www.tricemedical.com | customerservice@tricemedical.com